



EVO IPL, SHR, E-LIGHT MACHINE

- Hair removal
- Skin rejuvenation
- Pigment therapy
- Vascular therapy
- Acne treatment



THREE TREATMENT MODES

SHR IPL ELIGHT three treatment modes to meet a variety of treatment needs




INTELLIGENT OPERATING SYSTEM

The operator only needs to select the part and skin color to be treated, and the system will automatically generate the recommended treatment parameters. Even an inexperienced operator can easily grasp the operation of the instrument.





SHR "In-motion" Hair Removal

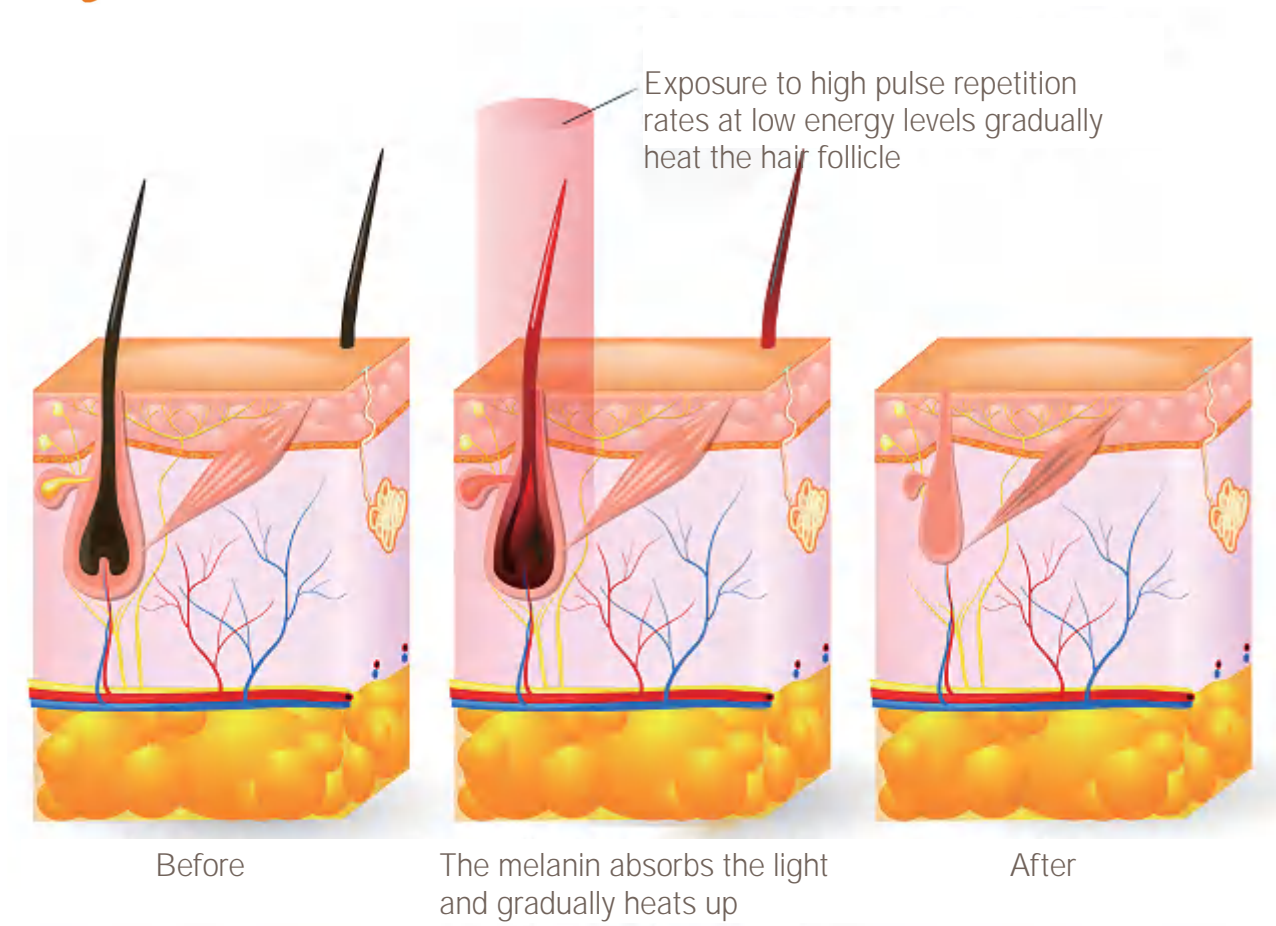

SHR stands for Super Hair Removal and it's the latest breakthrough innovation in IPL technology to launch onto the ever-advancing beauty and aesthetics market.



Faster

SHR in-motion treatment method has much improved the treatment speed. The SHR machine repetition rate is up to 10Hz. The handpiece could shoot 10 shots in 1 second.




More Comfortable

Unlike traditional IPL treatments, which use high energy levels that may cause discomfort and pain, SHR Skin Pulse treatments are mild and comfortable, SHR scans the treatment zone—gradually heating the skin with repetitive low fluence pulses, ensuring that you'll experience comfort even in the most sensitive of areas.




All Skin Types

The gradual heat distribution of SHR means treatments are totally comfortable, non-irritating and virtually pain free making it ideal for even sensitive dark skin types.

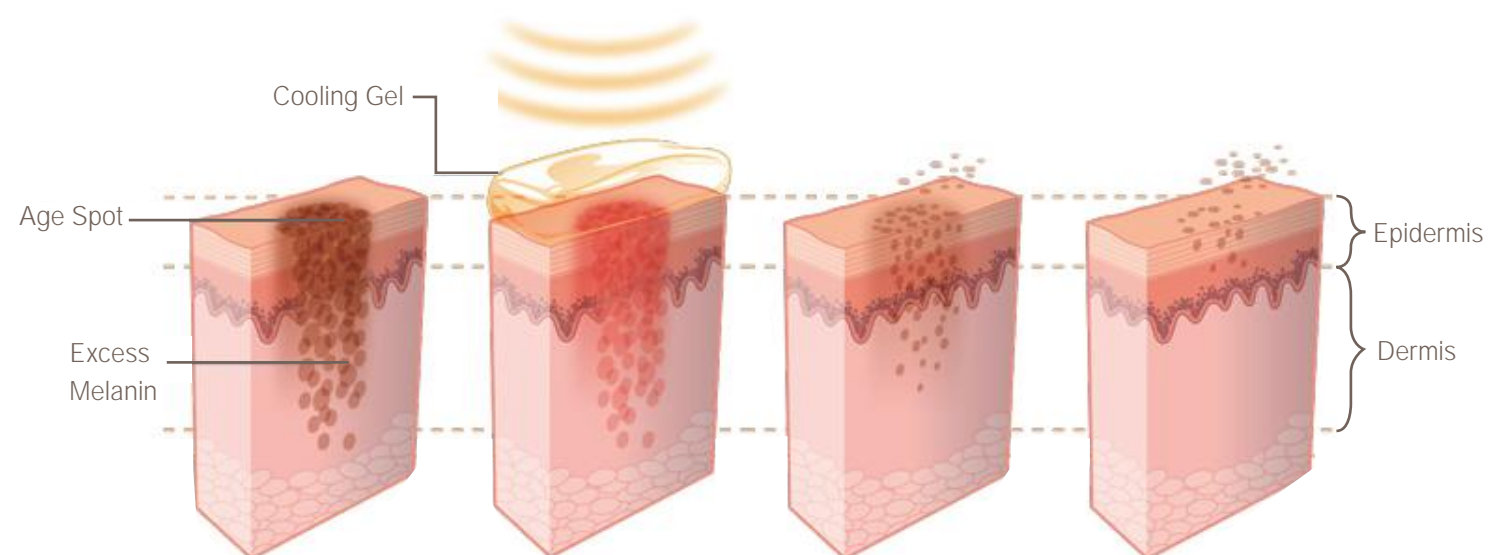

IPL technology is a type of Intense Pulsed Light that generates intense bursts of red or blue light. These heat the dermal tissue without harming or removing the outer epidermal layers of the skin. The body's natural processes then work to break down and remove these marks. The appearance of any skin imperfection is then visibly reduced.

Elight

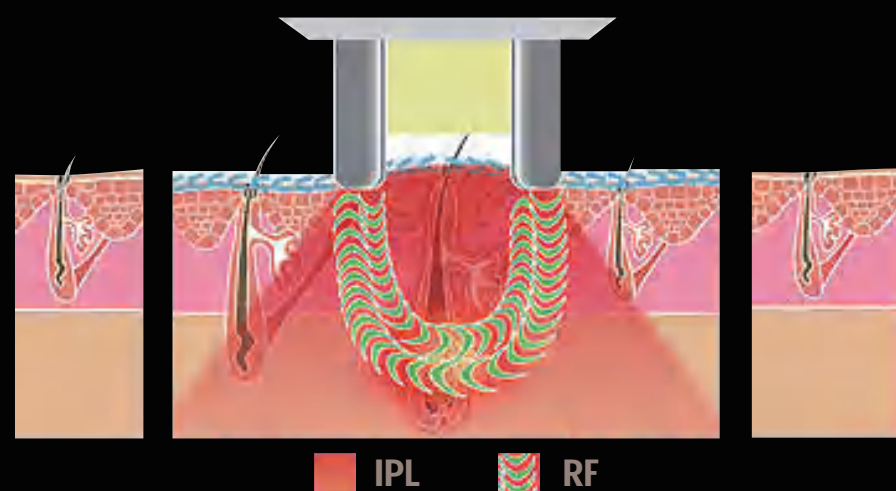

Elight is IPL+RF, while doing IPL treatment, radio frequency energy is output from the metal splints on both sides of the IPL treatment head, which has the effect of lifting and tightening





REAL-TIME INTELLIGENT MONITORING SYSTEM

20 intelligent detection functions, including water flow, water temperature, handle temperature, handle connection, water level etc.
Ensuring the stability of the machine and find out the possible problems of the machine in time



HR Handle

For in motion hair removal treatment with high power mode



Multi-function Handle

equipped with 6 different wavelength filters, which can be used for the treatment of acne, skin rejuvenation, spider veins, pigment etc.



handle connection status detection

Cooling sheet work detection, etc

Water level detection

2000W High power

The output energy is high enough to ensure that the energy will not decay even in the fast light mode, ensuring a good treatment effect.

handle connection status detection

Strong Cooling System


Ensuring handle surface is always cool and the treatment always comfortable.

At the same time, making sure that the machine works continuously for 24 hours without damage by the high temperature.

Water level detection

Six different filters

One handle is equipped with 6 filters with different wavelengths, which can meet various treatment needs

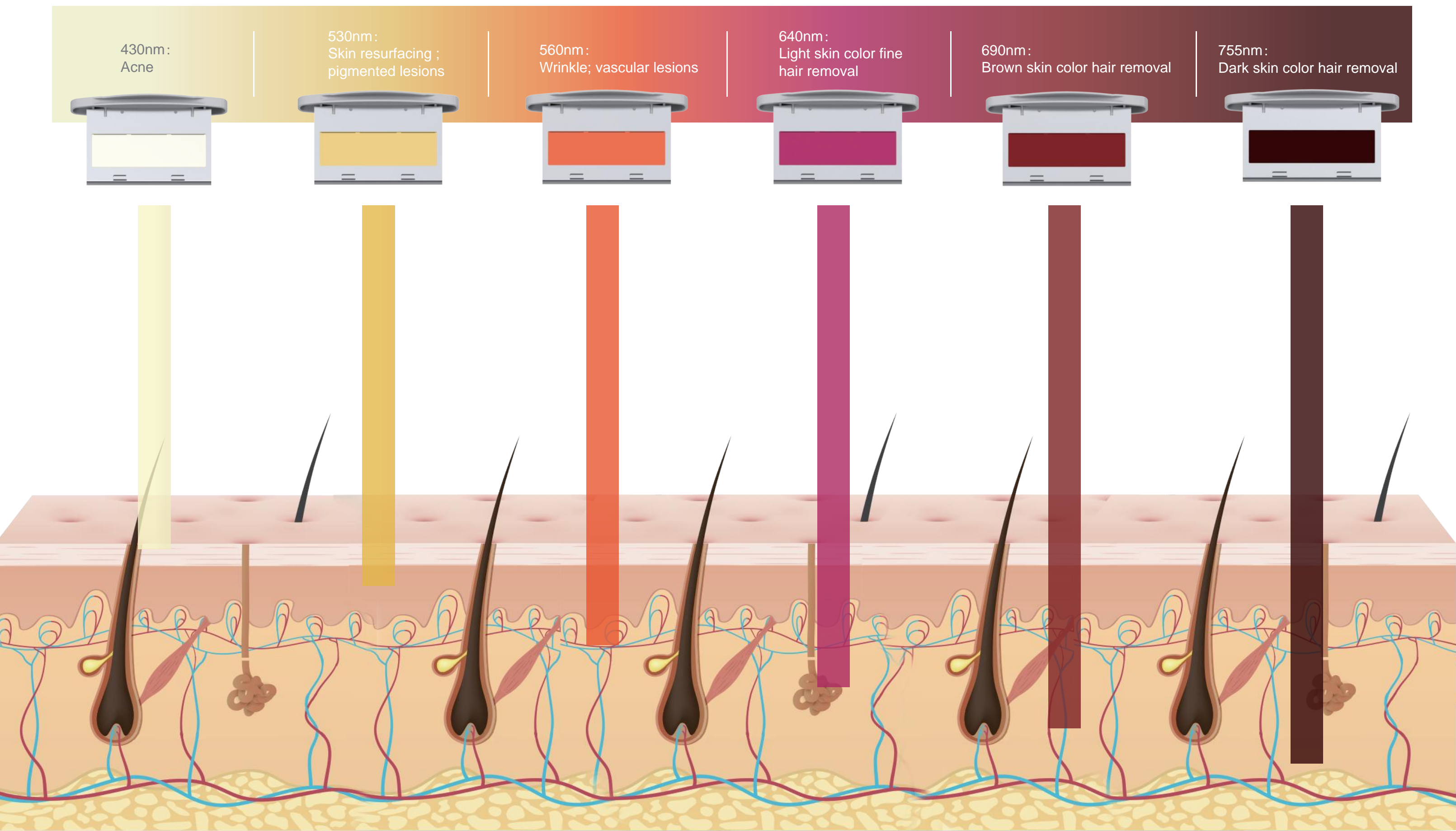



Handle Cooling

The handle system adopts 108W high-power refrigeration film, which is twice as much as other machines on the market, ensuring comfortable and safe treatment



Silver Reflector

The reflector is made of silver, which is not easy to be oxidized, has higher energy and longer life.


Sapphire Crystal

The handle crystal adopts sapphire crystal, which has better conduction of cold and can ensure the comfort of treatment



UK Xenon Lamp

The UK xenon lamp, diameter 9MM, high energy and long life. Warranty 1 million shots.

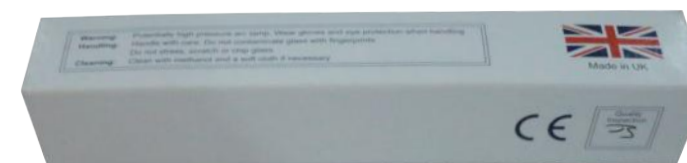



SPECIFICATIONS

Power	1500W
Wavelength	SR 560-950nm, HR 640-950nm
working modes	IPL, SHR, E-light
IPL energy	10-60j/cm2
Spot size: SR	8*40mm, HR: 15*50mm
SHR energy	3-15j/cm2
SHR frequency	1-10Hz
SHR pulse width	1-10ms
RF power	160W
RF energy	1-50j
Screen size	12.1 inches
Voltage	AC220V±10% 10A 50HZ , 110v±10% 10A 60HZ

